



Cap Illustrations

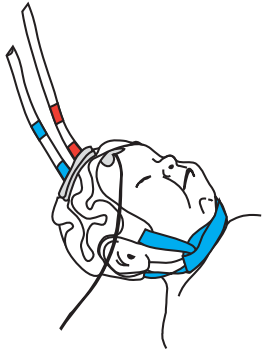


#1

Vinyl cooling cap is placed over scalp temperature sensor. Vinyl cap is secured under chin via Velcro fasteners.

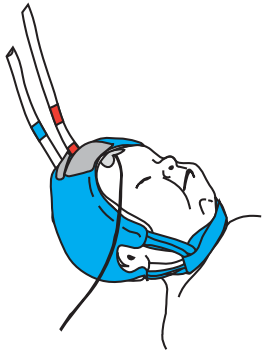


#2

Fabric retaining cap is placed over vinyl cap. Tubes route out top of fabric cap and are stabilized by Velcro closures. Fabric cap is secured under chin via vinyl cap's Velcro fasteners on left...

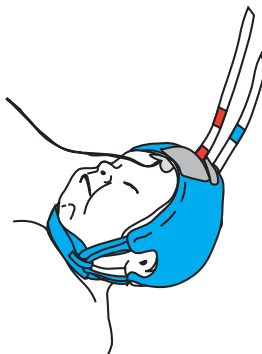

#3

...and right sides.

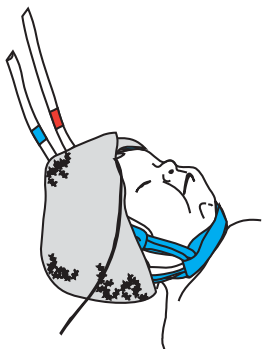

#4

Reflective silver thermal cap is placed over fabric cap. Silver cap attaches to fabric cap at three Velcro points (left, back, right). Tubes route out top of silver cap and are stabilized by Velcro closure.